



Pocket Guide

Babylog® VN800 / Babylog® VN600

Software version 1.n

Content

Product description	04
Operating concept	08
Getting started	12
Operation	16
Alarms	20
Monitoring	22

The Babylog® VN800 / VN600 Pocket Guide is not a replacement or substitute for the Instructions for Use. It is for informational purposes only. Any use of the device requires full understanding and strict observation of the Instruction for Use.

The availability of functionality depends on the country and availability. Options may not be included in a present device or may come at extra costs.


The pocket guide may not reflect the latest software version of the device.

Intended use

The intensive care ventilator Babylog is intended for the ventilation of neonates, infants, and children from 0.4 kg (0.88 lbs) up to 20 kg (44 lbs) or up to 6 years of age. This device provides mandatory ventilation modes and ventilation modes for supporting spontaneous breathing and also ventilation monitoring.

Product description


Intensive care Ventilator

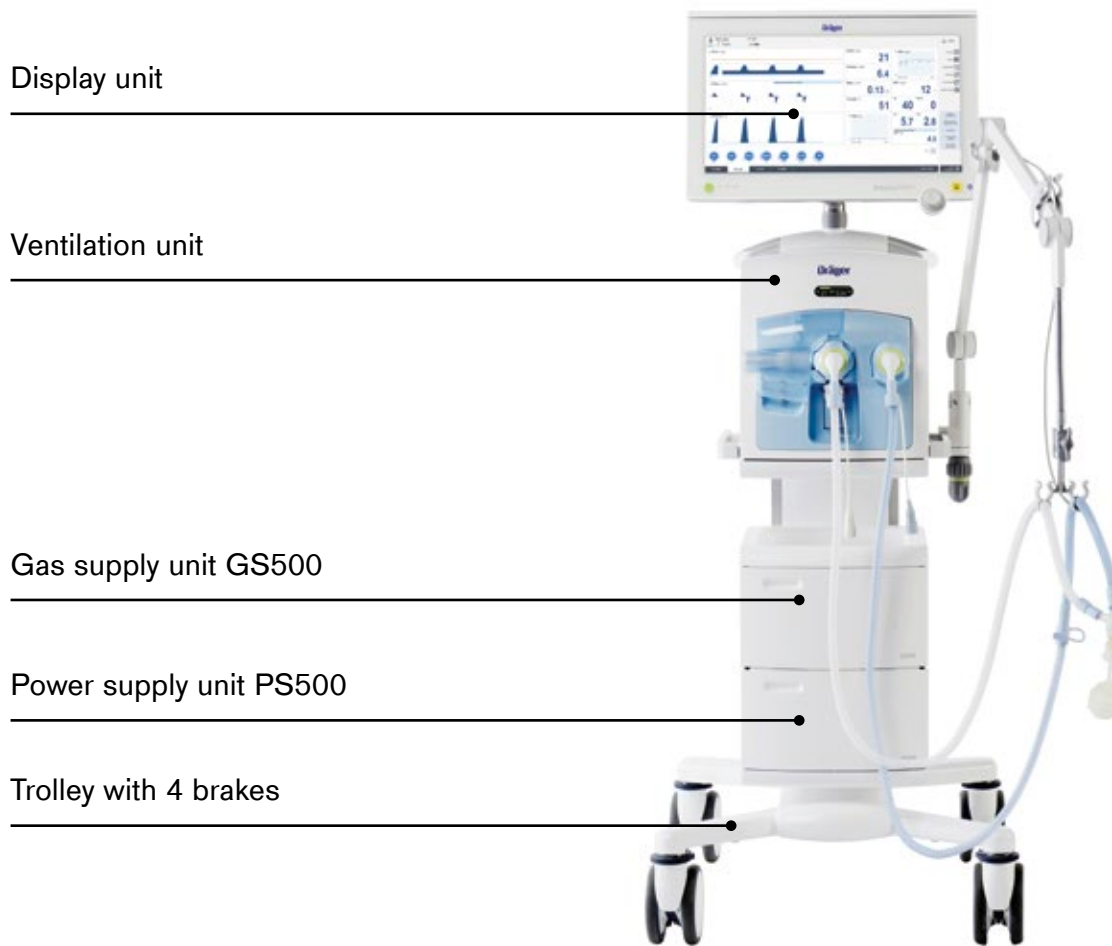

Babylog® VN800


Product description

Display unit, front

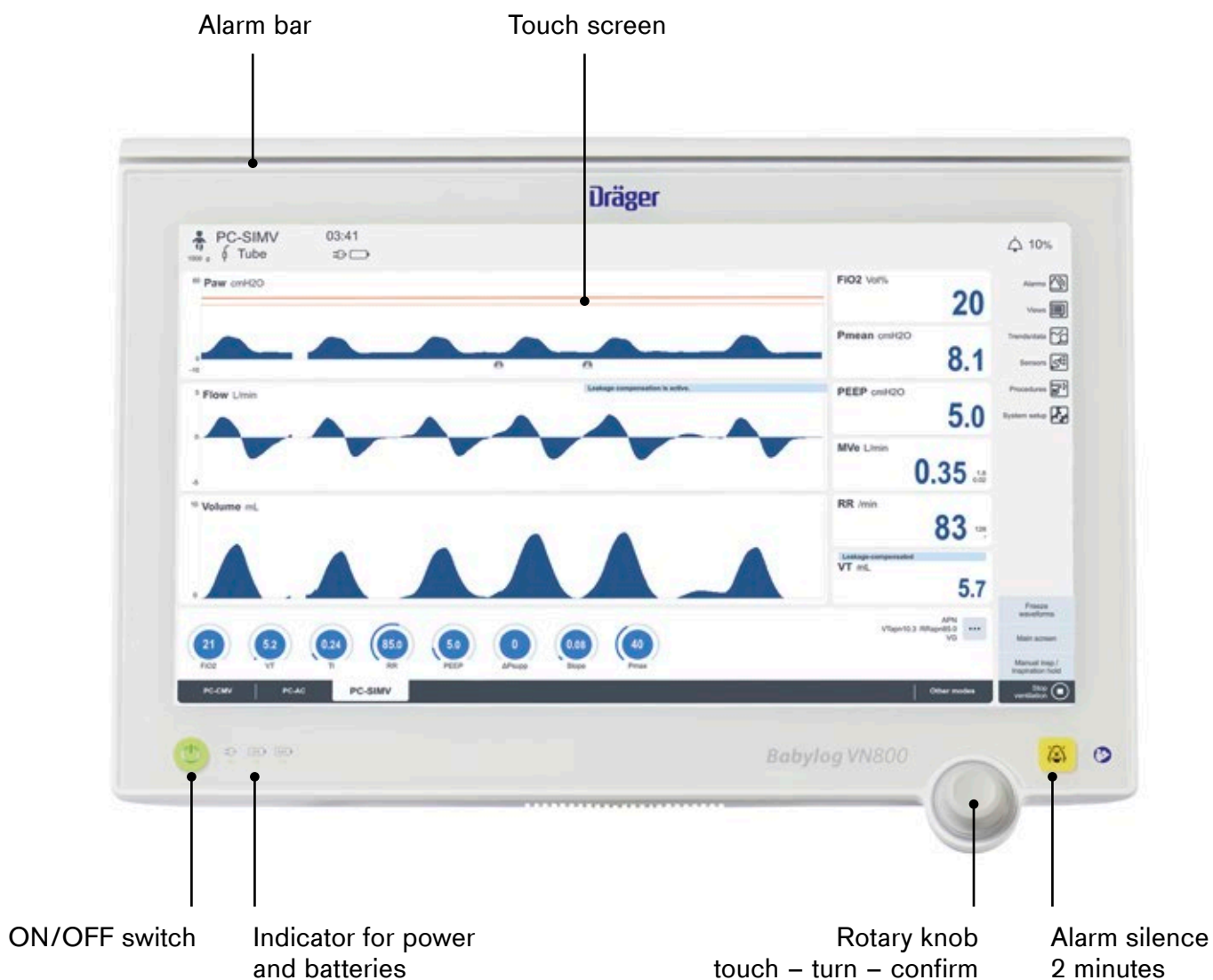


Product description

Display unit, rear

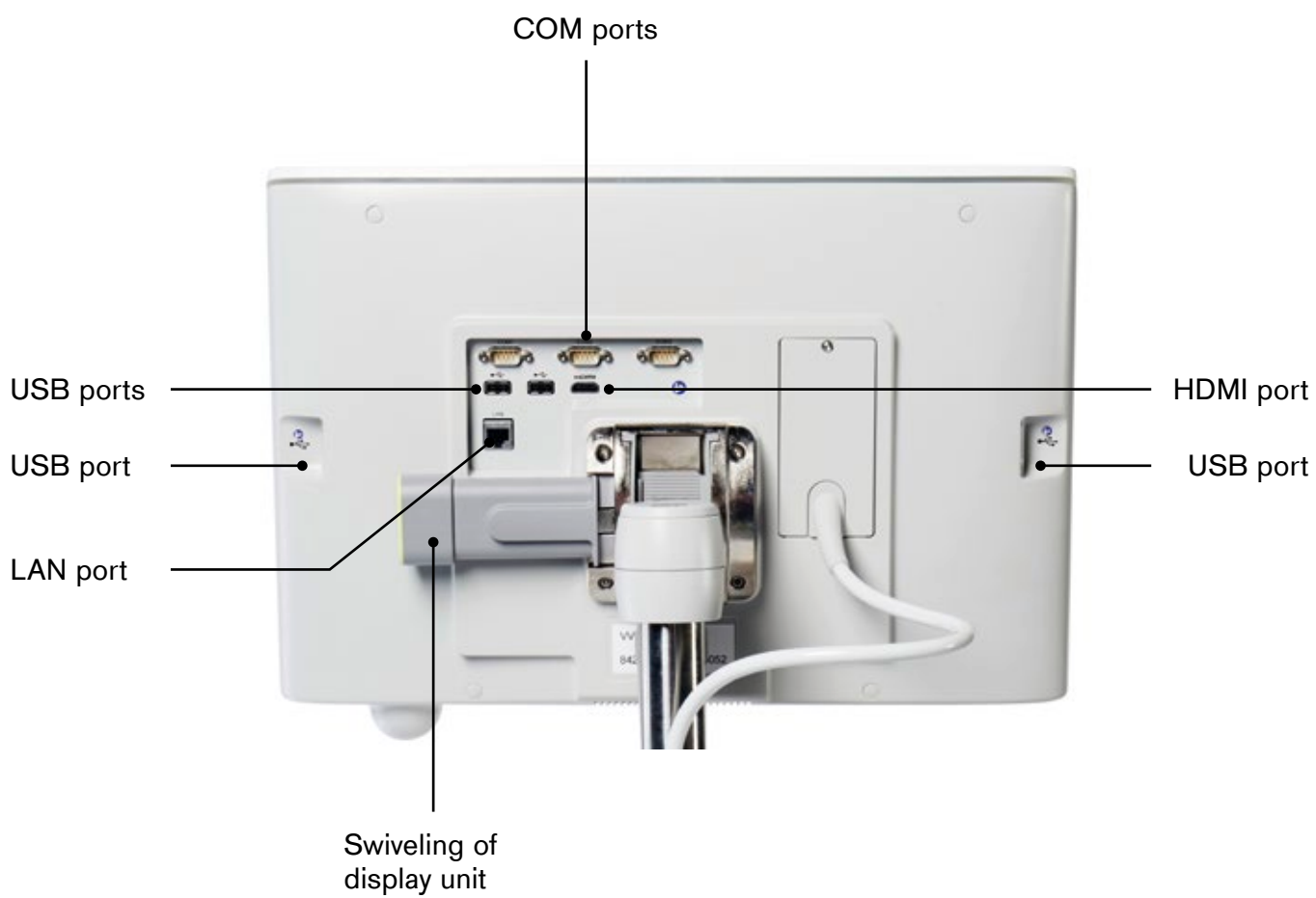

Product description

Ventilation unit, front

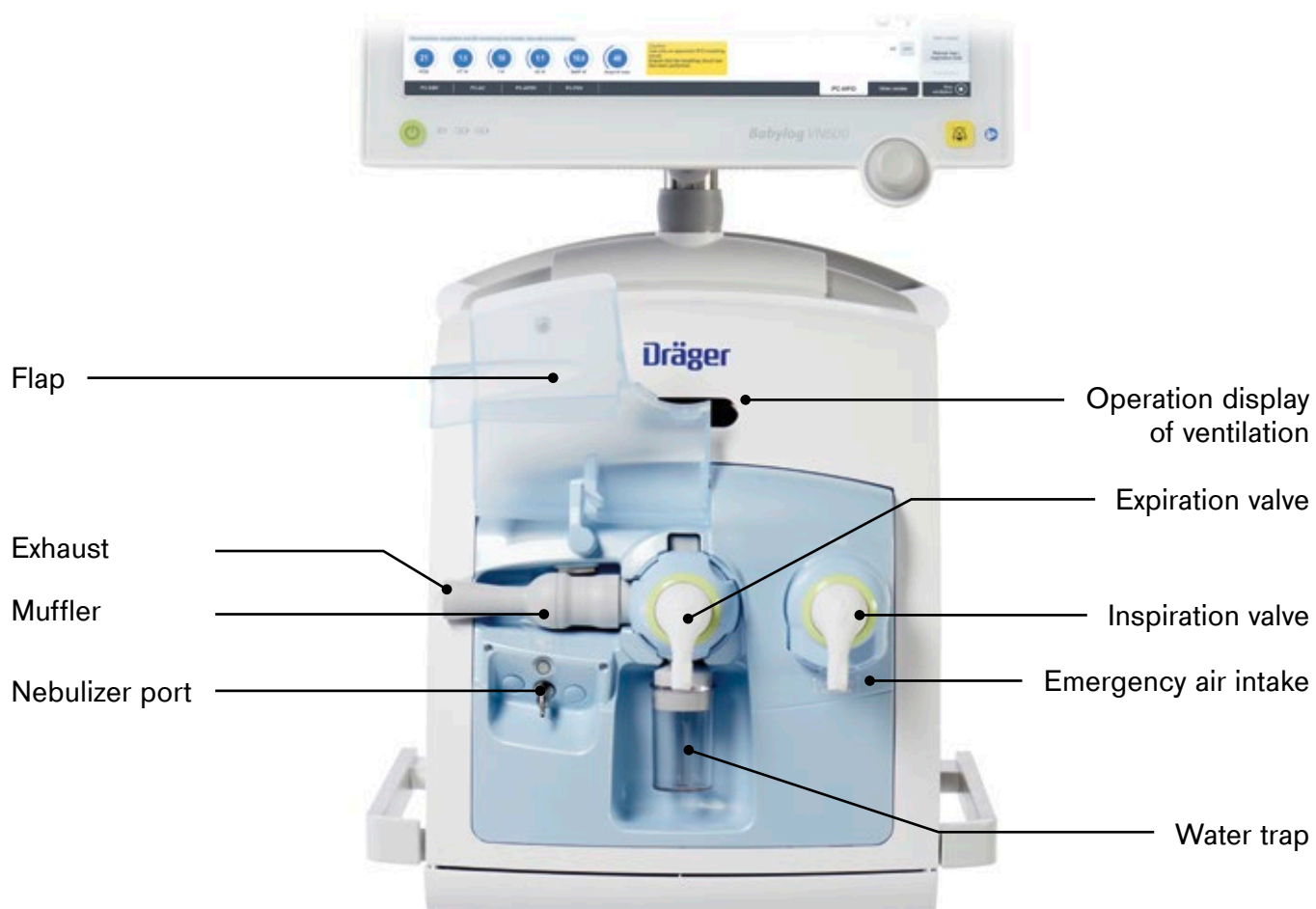


Operating concept

Screen, general structure



Operating concept

Main screen

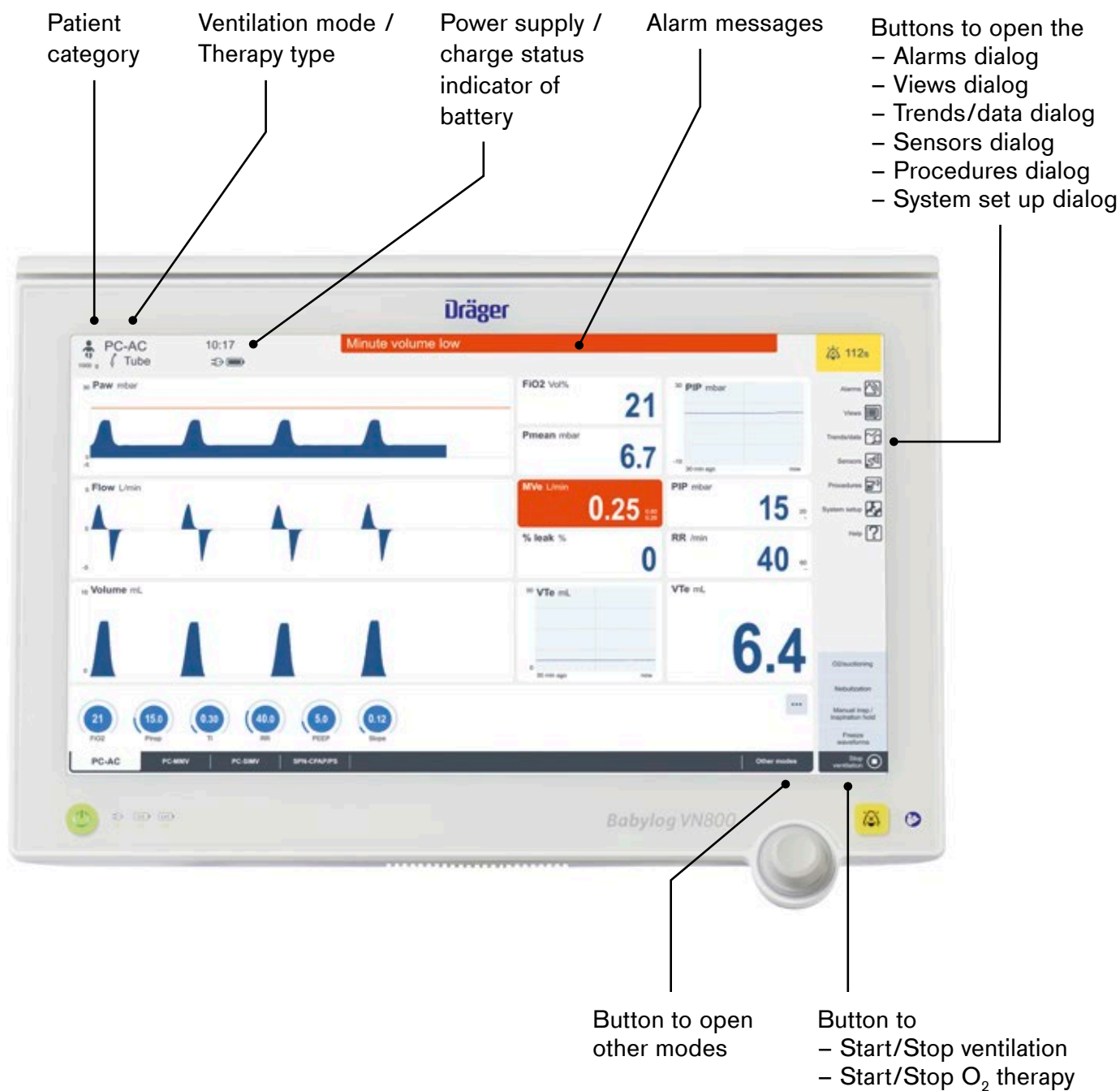

Operating concept

Using the screen – Colours

Alarms

The following colours are used to identify the priority of alarms:

Colour	Alarm priority	Action required
Red	High !!!	Immediate action is necessary in order to avert an acute danger
Yellow	Medium !!	Prompt action is necessary in order to avert a danger
Turquoise	Low !	Attention is necessary

Minute volume low



Tidal volume high



Battery activated




Operating concept

Using the screen – Colours

Rotary knob

If the device is turned on, the rotary knob is back lit with a colour.

- | | |
|--------|--|
| Orange | The selected function or setting must be confirmed by pressing the rotary knob. 10 seconds before the time during which the confirmation may be performed ends, the background lighting flashes. |
| Blue | After the device is turned on, the background lighting shines for 2 seconds. |

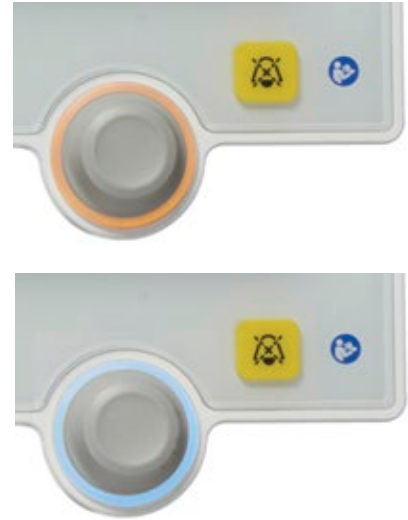

Touch – Turn – Confirm

Selecting and setting parameters and functions

Touch

the control element.
The control element turns orange.

Turn

the rotary knob to the right or to the left to make the setting.
The value is displayed in the control element.

Confirm

by pressing the rotary knob.
The setting is applied, and the control element turns gray or blue again.



Getting started

Turning on the device

Prerequisites:

- The accessories needed for operations on the patient are connected.
- The mains power supply and the gas supply are connected.
- The main switch is turned on.

Ventilation unit,
left side:
Main switch



1. Press the **ON/OFF** key on the display unit. The system is started. The start dialog is displayed.
2. Touch the **New patient** button if a new patient is admitted.
or
3. Touch the **Current patient** button to continue with existing settings.



Getting started

System test

Starting the system test

1. Touch the **Start** button in the System test line.
2. Confirm with the **rotary knob**.
3. Answer the questions in the information field by touching the **Yes** or **No** button.


After a successful system test has been performed, a message appears for continuing the **breathing circuit test**. Proceed with 'Yes'.



Breathing circuit test

Performing the Breathing circuit test


1. Touch the button in the Breathing circuit test line.
2. Select the **breathing circuit** used from the pick list and the type of **humidification**.
3. Confirm with the **rotary knob**.
4. Follow the instructions in the information field and touch the button **OK** when asked to do so.

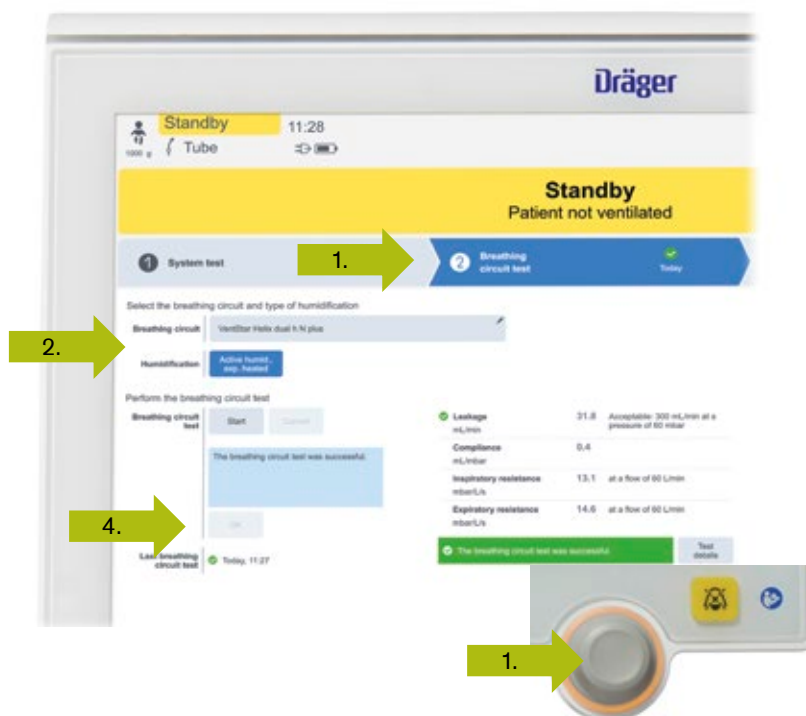

Getting started


Starting the therapy (based on the patient height)

1. Touch the button for the relevant **patient category**.
2. Confirm with the **rotary knob**.
3. Touch the button for **weight**.
4. Set the body height by turning the **rotary knob** and push to confirm.
5. Touch the button for the relevant **therapy type (Tube, NIV, O₂ therapy)**.
6. Confirm with the **rotary knob**.
7. Check the alarm limits and touch the button **Start ventilation**.
8. Confirm with the **rotary knob**.

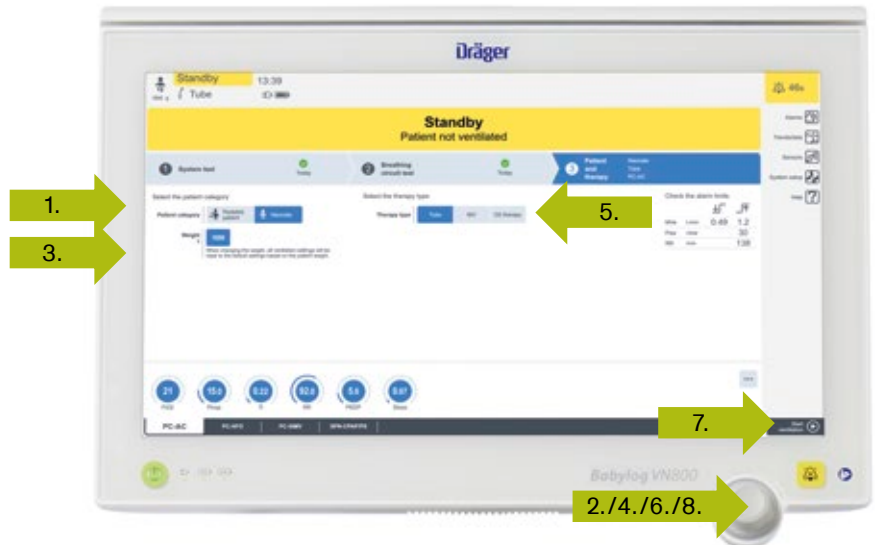

Starting the therapy (based on the patient category)

1. Touch the button for the relevant **patient category**.
2. Touch the button for the relevant **therapy type (Tube, NIV, O₂ therapy)**.
3. Confirm with the **rotary knob**.
4. Check the therapy settings (ventilation mode, ventilation parameter, alarm limits) and touch the button **Start ventilation**.
5. Confirm with the **rotary knob**.

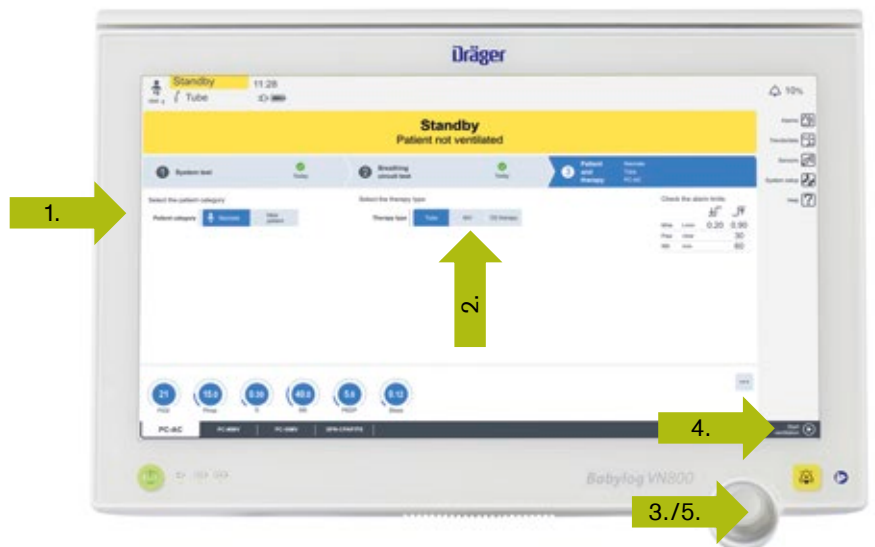


Getting started

Therapy Quick start ventilation

In an emergency situation, therapy can be started immediately without operational readiness being checked.

1. Touch the button **Quick start ventilation**. The patient and therapy step is displayed.
2. Press the **rotary knob**. The therapy is started. The main screen is displayed.

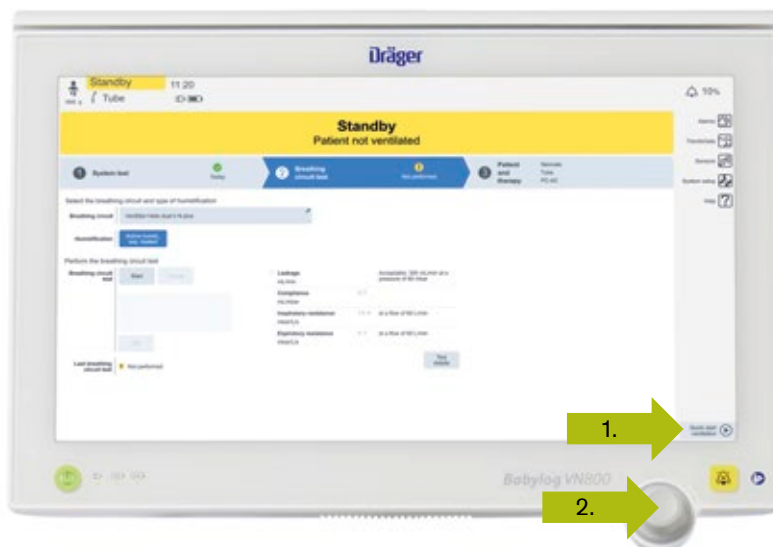

Operation

Changing the ventilation mode

1. Touch the relevant tab.
2. The colour of the tab turns orange (if necessary, adjust the ventilation parameter).
3. Confirm with the **rotary knob**.
The color of the tab turns white.



Setting ventilation parameters

1. Touch the corresponding **therapy control** (colour turns orange).
2. Set the value by turning the **rotary knob** and push to confirm.



Operation


Advanced settings

1. Touch the button **Advanced settings**. The advanced settings of the selected ventilation mode are displayed and can be configured.
2. The current settings are displayed in an information field next to the button.



Changing the therapy type during therapy

1. Touch the field for the therapy type in the header bar.
2. Touch the button of the relevant therapy type
 - Tube
 - NIV
 - O₂ therapy
3. Confirm with the **rotary knob**. The device is in standby mode.
4. Start the therapy by touching the **Start ventilation** button.

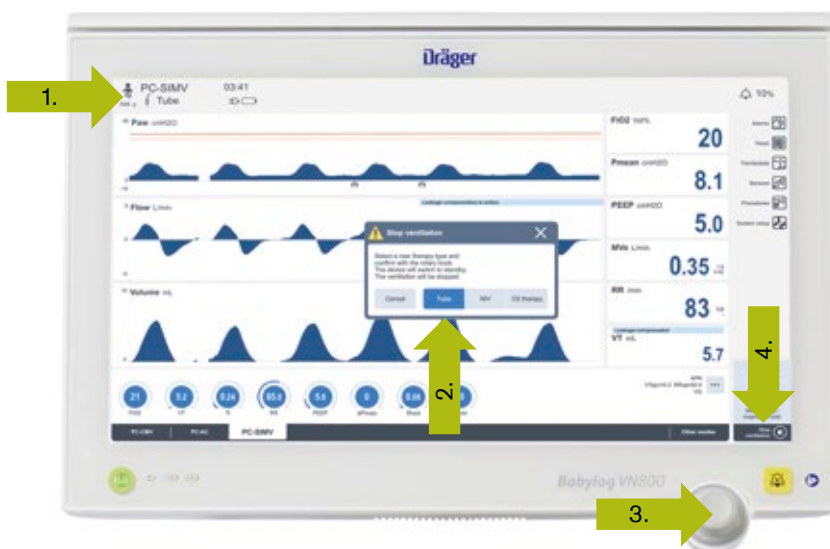

Operation

Non-invasive ventilation (NIV)

In the **Adult** and **Pediatric patient** categories, all the ventilation modes can be selected. In the **Neonate** patient category, only the **SPN-CPAP** and **PC-CMV** ventilation modes can be selected.



O₂ therapy

Patient-side accessories such as oxygen masks or nasal cannulas can be used for O₂ therapy.

During O₂ therapy, only the O₂ concentration, the inspiratory flow, and the inspiratory pressure are monitored.

The alarm limits for the parameters **VT**, **MVe**, **RR**, **Paw**, and **Tapn** are not active. The alarm limits for O₂ monitoring are automatically set by the device.



Operation

Interrupting therapy

1. Touch the button **Stop ventilation** in the therapy bar. A pop-up window opens.
2. Touch the button **Stop ventilation**.
3. Confirm with the **rotary knob**.
4. The following message is displayed in the header bar: Standby mode activated.
5. Touch the **Alarm reset** button in the header bar.
6. Confirm with the **rotary knob**.




Alarms

Display of alarms – Optical and acoustic alarm signals



Display of alarms – Confirming alarm messages

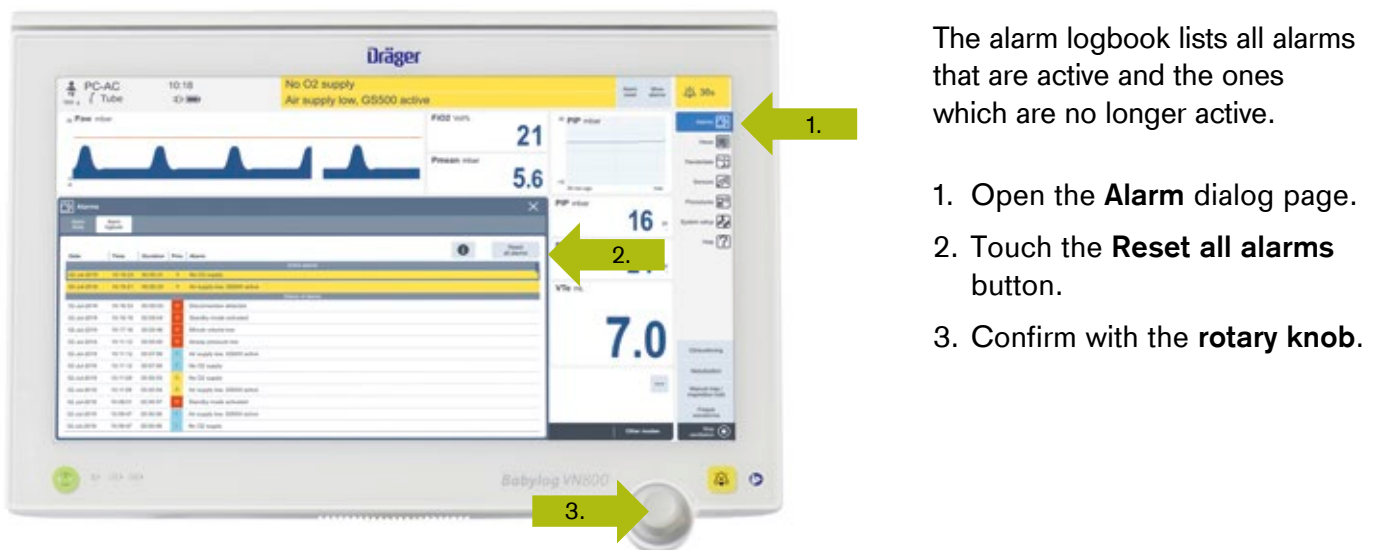

Alarms

Set alarm limits

1. Touch the **Alarms** button in the main menu bar.
The **Alarm limits** dialog page is preselected.
2. Touch the corresponding button for the upper or lower alarm limit.
3. Set the value by turning the **rotary knob** and confirm.



Neonatal flow monitoring

The flow sensor must be manually calibrated at least once a day and in the following cases as well:

- During the system test and before use
- After replacing the following components:
 - Sensor insert
 - Flow sensor

If the flow sensor was unplugged only briefly, it does not need to be recalibrated.



1. Port for the neonatal flow sensor




Neonatal flow sensor

Neonatal flow monitoring

Calibrating the neonatal flow sensor

Performing the calibration, open the dialog page

1. Touch the **Sensors** button in the main menu bar.
2. Touch the **Neonatal flow sensor** tab.

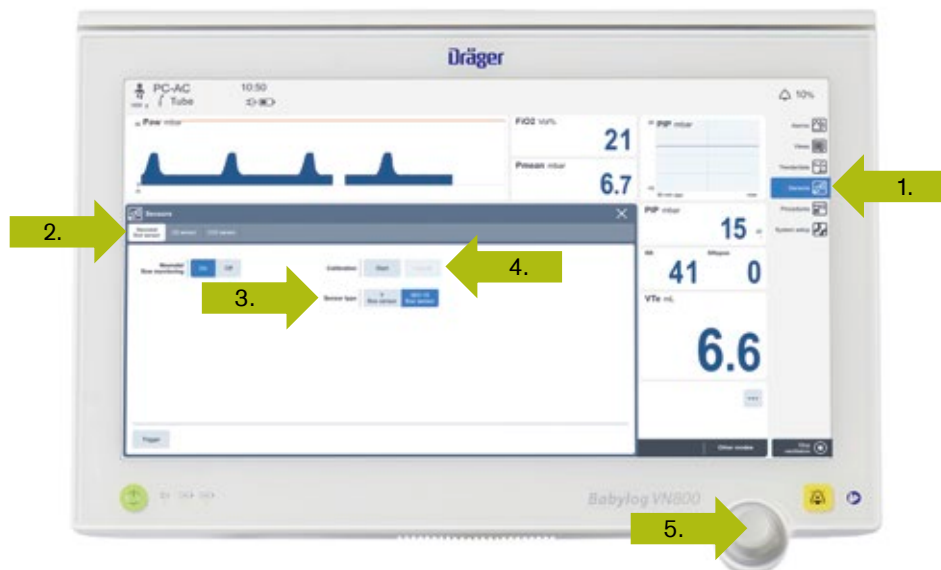

3. Select the sensor type used in the line **Sensor type**.

4. Touch the **Start** button.

Perform the instructions displayed in the tooltip.

5. Confirm the instructions with the **rotary knob**.

The following neonatal flow sensors are available:



Flow sensor Y-piece



Flow sensor ISO 15

Perform the instructions displayed in the tooltip:


- Remove the tube connector
- Seal the flow sensor

5. Start the calibration by pressing the **rotary knob**.

The tooltip displays information.

After the calibration is complete, the **Start** button turns blue.

After calibration of the neonatal flow sensor connect the tube connector.

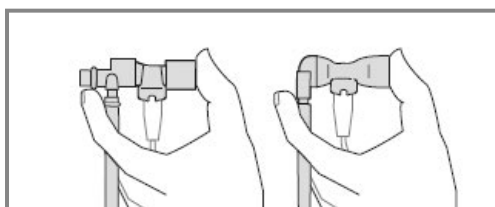

Sealing the neonatal flow sensor
(wear a sterile glove when sealing the flow sensor)

Not all products, features, or services are for sale in all countries.
Mentioned Trademarks are only registered in certain countries and not necessarily in the country in which this material is released. Go to www.draeger.com/trademarks to find the current status.

CORPORATE HEADQUARTERS

Drägerwerk AG & Co. KGaA
Moislinger Allee 53–55
23558 Lübeck, Germany

www.draeger.com

Manufacturer:

Drägerwerk AG & Co. KGaA
Moislinger Allee 53–55
23542 Lübeck, Germany

Locate your Regional Sales
Representative at:
www.draeger.com/contact



REGION EUROPE

Drägerwerk AG & Co. KGaA
Moislinger Allee 53–55
23558 Lübeck, Germany
Tel +49 451 882 0
Fax +49 451 882 2080
info@draeger.com

REGION MIDDLE EAST, AFRICA

Drägerwerk AG & Co. KGaA
Branch Office
P.O. Box 505108
Dubai, United Arab Emirates
Tel +971 4 4294 600
Fax +971 4 4294 699
contactuae@draeger.com

REGION ASIA PACIFIC

REGION ASIA PACIFIC
Draeger Singapore Pte. Ltd.
61 Science Park Road
The Galen #04-01
Singapore 117525
Tel +65 6872 9288
Fax +65 6259 0398
asia.pacific@draeger.com

REGION NORTH AMERICA

Draeger, Inc.
3135 Quarry Road
Telford, PA 18969-1042, USA
Tel +1 800 4DRAGER
(+1 800 437 2437)
Fax +1 215 723 5935
info.usa@draeger.com

REGION CENTRAL

AND SOUTH AMERICA

Dräger Indústria e Comércio Ltda.
Al. Pucurui - 51 - Tamboré
06406-100 - Barueri - SP
Tel. +55 (11) 4689-4900
relacionamento@draeger.com