



DID YOU KNOW – LEAK COMPENSATION



Dräger

LEAK COMPENSATION

What is it?


Leak compensation is an advanced algorithm to precisely calculate the tidal volume that flows to the patient.

The goal is to display the leakage-corrected values of flow and volume waveforms and tidal and minute volume values.


Leak compensation of flow and volume curves and the volume measurements can be switched ON or OFF in the system setup.

How is it working?

If Leak Compensation is active the ventilator calculates a lung tidal volume (labelled as VT) using a sophisticated algorithm. The calculated lung tidal volume is not measured, but calculated, taking into consideration leaks.

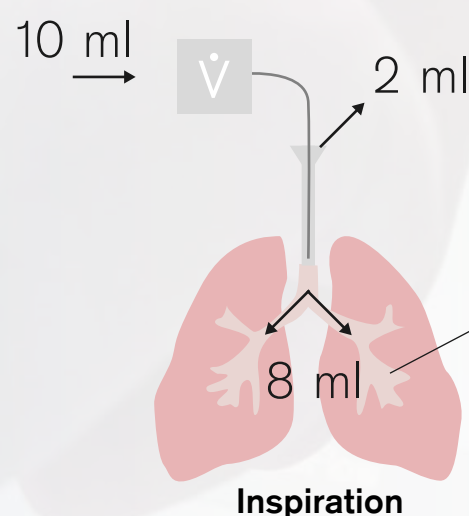

Leak Compensated VT
= VT in the lung

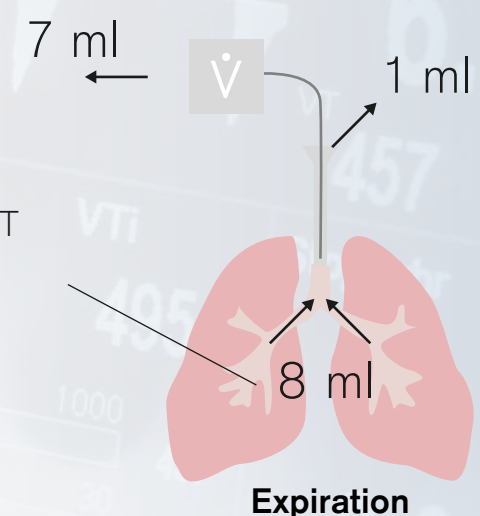

Why is it helpful to improve outcome?

Leak compensation (even together with leak adaptation) is a powerful combination to **decrease the users' cognitive workload and focus on the patient treatment and outcome.**

Leak compensation might reduce VT low alarms.

The user can monitor the inspiratory and expiratory tidal volume and the calculated tidal volume in the lung.

VT Alarms are based off the VT instead of VT_e which **may decrease the number of alarms during high leak situations** but it's important to monitor the VT_i, VT_e and the VT to ensure effective ventilation is delivered to the patient.

