

THIS IS VENTILATION

NEONATAL VENTILATORS

Babylog® VN800

Babylog® VN600



Our mission: Protecting fragile lungs & the developing brain



14-38%

preventable deaths due to Broncho-pulmonary Dysplasia in ELBW infants*.

We support you to provide protective care and an environment that fits the baby's needs. We want to help you to reduce the number of respiratory complications and preventable deaths so as to improve your patient outcome.

OPTIMISING OXYGENATION AND VENTILATION

The goal is to keep tidal and minute volume stable for an appropriate gas exchange and to enable lung and brain protective ventilation. In order to optimise your patient's oxygenation, strategies and ventilation modes are needed to prevent fluctuations in carbon dioxide levels as well as stabilising blood gases. A developmental care friendly environment helps to improve the outcome of neonates as well.

ENABLING QUICK & SAFE WEANING

The start of ventilation therapy is always the start of a weaning process. To prevent chronic lung disease such as bronchopulmonary dysplasia (BPD), it is essential to directly help the infant regain their own respiratory drive. We support the earliest possible weaning from mechanical ventilation. With all of our ventilation applications we want to reduce the work of breathing and encourage spontaneous breathing.

SUPPORTING INTUITIVE OPERATION

"Know your tool" – in order to enable lung and brain protective ventilation, you can use different modes – but these modes can't be used to their full potential unless you are able to monitor, analyse and adapt parameters. Configurable and meaningful views help with intuitive operation which lowers education times and helps to avoid human errors.

* Davidson LM, Berkelhamer SK: Bronchopulmonary Dysplasia: Chronic Lung Disease of Infancy and Long-Term Pulmonary Outcomes. J Clin Med 2017, 6(1):4. doi: 10.3390/jcm6010004.



All ventilation strategies combined in one device




Babylog® VN800



Babylog® VN600

HIGH-END NEONATAL VENTILATION

The Dräger Babylog comes with a new user interface and design which makes it easier and safer for you to operate. The neonatal ventilator supports lung and brain protective ventilation modes throughout the whole respiratory cycle and can be easily integrated in a developmental care-friendly workplace.

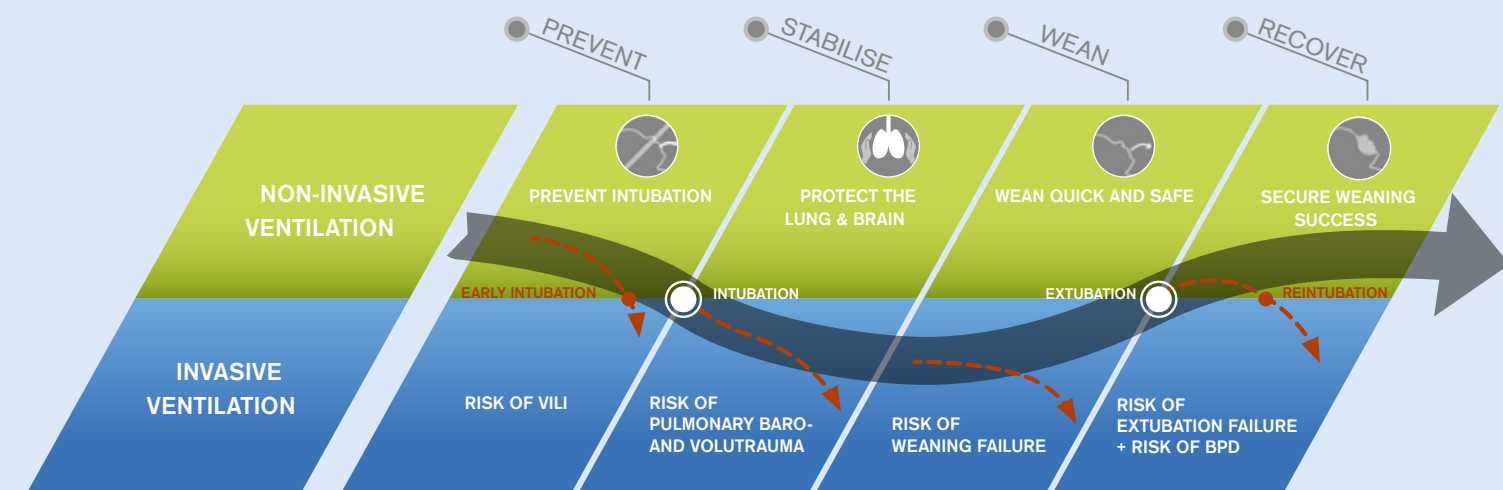

RESPIRATION PATHWAY

Neonatal ventilation is an important part of the care delivered in the neonatal intensive care unit. The aim of any ventilation strategy is to support the premature infant's respiratory system without inducing damage to the lung or the brain. Those ventilation strategies can be viewed as a continuum of dependencies that accompany the neonate and the caregiver from beginning to the end of respiratory challenges as pictured in our respiration pathway.

PREVENT

Avoiding intubation helps to protect newborn's lungs and brain. Only when necessary, should it be done.

STABILISE

Strategies that include synchronised respiratory support and volume guarantee ventilation allow you to stabilise your ventilated infant.

WEAN

Dedicated therapy approaches facilitate and automatise weaning.

RECOVER

With non-invasive ventilation, you minimise work of breathing as well as improve the healthy recovery in ventilated infants.

Learn more about the respiration pathway and understand the different strategies and treatment objectives for neonatal ventilation: draeger.com/neonatal-ventilation

Clear user interface

that guides you in a hectic environment




OPERATING PRINCIPLE

- Quick and safe to operate even in the most stressful situations due to intuitive menu access to both settings and your clinical data.
- All patient data, alarms and trends are fully recorded. Conveniently exported via USB interface.
- Switch between multiple view configurations with the touch of a finger.
- Step-by-step guidance leads you through every procedure.


A BRILLIANT USER INTERFACE

- Easy to read and navigate thanks to our new colour concept and glass touch display.
- Safe and intuitive operation is supported by tablet-like operating system and easy access to sub menus.
- The 360° alarm light flashes in the colour of the corresponding alarm priority and is visible from every direction.

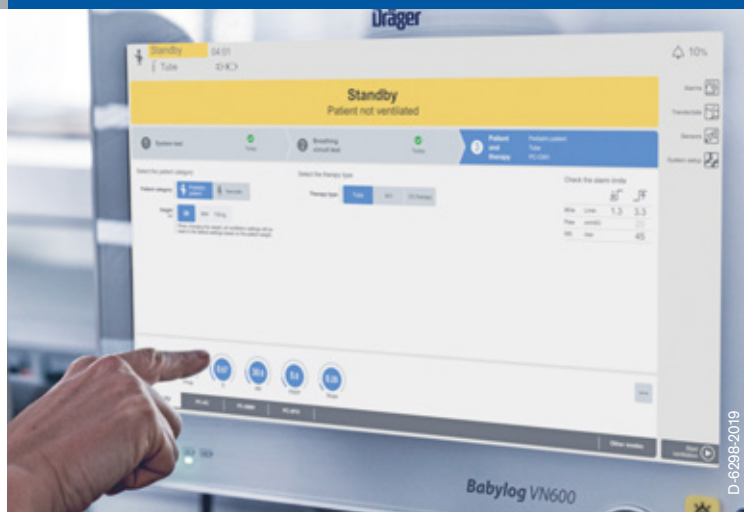

PATIENT TRANSPORT

- No need to change the ventilator during patient transport.
- Operate the device in intra-hospital transport independently with gas and external power supply units.
- Utilise the bed coupling to safely transport the patient in the incubator or bed.

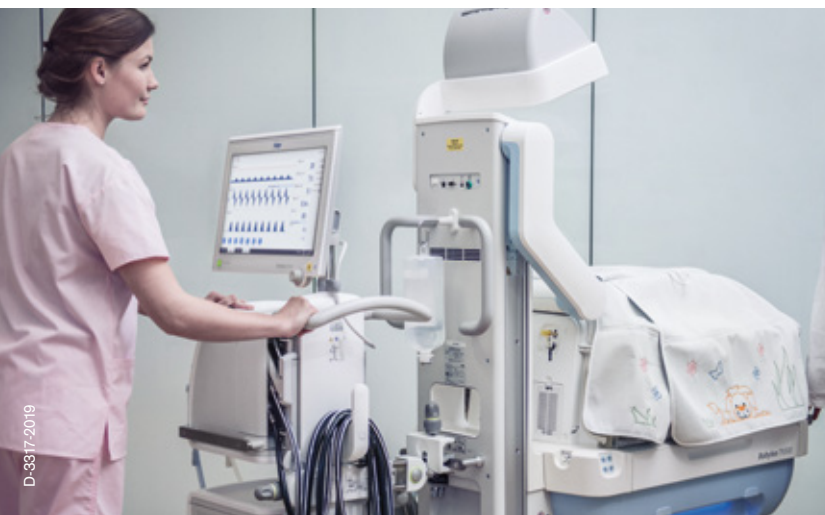

Lung and brain protective care



LUNG AND BRAIN PROTECTIVE VENTILATION

Our set of treatment tools support you in applying the right protective ventilation strategy in order to prevent lung injury, haemodynamic and neurological impairment.

- Dedicated invasive and non-invasive ventilation capabilities including high-flow oxygen therapy.
- Lung and brain protective ventilation due to automated pressure regulation with the original Dräger Volume Guarantee.
- Lung and brain protective ventilation due to High Frequency Ventilation with Volume Guarantee (HFO-VG).
- Stable minute ventilation and protective weaning with Mandatory Minute Ventilation (PC-MMV/VG+PS).
- Maintain reliable and sensitive triggering and stable lung tidal volumes with original Dräger leak adaptation and leak compensation technology.
- Proportional support for compensation of ETT resistances.

FURTHER PROTECTIVE THERAPY SOLUTIONS:

Thermoregulation

Avoiding cold or heat stress is one of the major targets for thermoregulation applications on the NICU. Actively providing a stable climate for premature infants, helps them to thrive. Dräger supports you in improving your clinical outcomes with open care warmers, closed care incubators, and combined solutions.

Jaundice Management

The avoidance of invasive interventions reduces stress and pain for the baby. Our non-invasive screening and therapy devices allow for effective jaundice management to help you minimise jaundice risks at an early stage while also reducing stress for the newborn.

Developmental care friendly environment

CARE-CENTRED WORKPLACES

From delivery to discharge: As your specialist in acute care, we want to accompany you through your clinical patient pathway and enable a developmental care- friendly environment to support all of the complex needs of the developing lung, brain and other organs. We do this through our wide range of products and solutions for L&D, transportation, and the NICU. Our products:

- Are compatible with each other and work with the same Dräger operating philosophy.
- Give you a flexible workplace integration of Babylog with different cockpit sizes and mounting possibilities.
- Offer low operational noise levels even during High Frequency Ventilation or non-invasive ventilation.
- Provide you with effective infection prevention thanks to easy cleaning of the glass touch screen and other smooth surfaces.
- Come with longer circuits and cables which allow staff and parents to remove the baby from the incubator for kangaroo time without compromising the baby's ventilation.
- Include a broad range of Dräger accessories with optimised circuits for High Frequency Ventilation and the non-invasive interface Babyflow plus.



10%
of all babies worldwide are born
prematurely and need dedicated care.*

* Website WHO organization



Connectivity and comprehensive services

CONNECTIVITY

We envision a future of acute care where medical devices are connected as a system. Interoperability between different devices can help to avoid preventable medical errors and potentially serious inefficiencies.

The new standardised network protocol named SDC makes the safe and dynamic connectivity in the hospital possible which will allow the interoperability of medical devices in the future.

Our first step is connectivity through CC300:

- Full HL7 data export to HIS from all devices: reliably exchange high-quality data in a standardised format between medical devices and EMR.
- Future-proof open connectivity: standardised and secure communication between medical devices with a high level of cyber security.



D-2662-2018

360° SERVICES FOR THE COMPLETE PRODUCT LIFECYCLE

To ensure the high performance of our ventilators along the entire lifecycle, we provide exclusive knowledge and skills when installing and maintaining Dräger systems tailored to the specific needs of your hospital. We offer service that goes way beyond classic device maintenance. It encompasses comprehensive services prior, during and after the installation of your purchased units:

- Product Service: e.g. device maintenance
- Professional Service: e.g. IT consulting and system integration
- Training on our products and services: e.g. application training
- Multivendor Service: maintenance for your entire medical equipment regardless of the manufacturer
- Digital Services: e.g. network-based services and analysis of device data (not all service products are available in every country)

REFERENCE SITE

“Babylog’s new user interface is clear, concise and as you know it – well-structured”

Dr. Anja Stein, Senior Physician of Neonatology at the University Hospital Essen



D-5883-2019

TRAINING

To facilitate protective ventilation strategies in neonates, we provide practical tips and hands-on strategies to help you set-up, maintain and optimise your technology.



D-3246-2019

ACCESSORIES

Dräger original neonatal ventilation accessories are suitable for different ventilation strategies and have an optimised length, pressure and support humidification. Our accessories are completely compatible with our neonatal ventilator Babylog. Find them all in our neonatal accessories catalogue.



D-15747-2018

BABYFIRST.COM

Clinical knowledge and family focus – all in one place. BabyFirst is a clinical and informational knowledge base for neonatal professionals and parents of premature babies. It provides pertinent neonatal clinical information on various topics and it is a valuable tool for parents that are thrust into what can be a very emotional and unfamiliar situation and are looking for more information. Visit www.babyfirst.com

For more information, visit our website draeger.com/neonatal-ventilation



D-15337-2017



D-5561-2019

Not all products, features, or services are for sale in all countries.

Mentioned Trademarks are only registered in certain countries and not necessarily in the country in which this material is released. Go to www.draeger.com/trademarks to find the current status.

CORPORATE HEADQUARTERS

Drägerwerk AG & Co. KGaA
Moislinger Allee 53–55
23558 Lübeck, Germany

www.draeger.com

Manufacturer:

Drägerwerk AG & Co. KGaA
Moislinger Allee 53–55
23542 Lübeck, Germany

Locate your Regional Sales
Representative at:
www.draeger.com/contact



REGION DACH

Drägerwerk AG & Co. KGaA
Moislinger Allee 53–55
23558 Lübeck, Germany
Tel +49 451 882 0
Fax +49 451 882 2080
info@draeger.com

REGION EUROPE

Drägerwerk AG & Co. KGaA
Moislinger Allee 53–55
23558 Lübeck, Germany
Tel +49 451 882 0
Fax +49 451 882 2080
info@draeger.com

REGION MIDDLE EAST, AFRICA

Drägerwerk AG & Co. KGaA
Branch Office
P.O. Box 505108
Dubai, United Arab Emirates
Tel +971 4 4294 600
Fax +971 4 4294 699
contactuae@draeger.com

REGION ASIA PACIFIC

Draeger Singapore Pte. Ltd.
61 Science Park Road
The Galen #04-01
Singapore 117525
Tel +65 6872 9288
Fax +65 6259 0398
asia.pacific@draeger.com

REGION CENTRAL

AND SOUTH AMERICA

Draeger Panamá Comercial
S. de R.L.
59 East Street, Nuevo Paitilla,
House 31, San Francisco Town
Panama City, Panama
Tel +507 377-9100 / Fax -9130
servicioalcliente@draeger.com