



Dräger Atlan A300/A350 Ceiling and Wall Anaesthesia Workstation

The Atlan family's platform design offers floor-free anaesthesia workstations with ceiling and wall-mounted versions. They support organised clinical environments by enhancing workflow efficiency, improving patient access, simplifying cleaning procedures, and optimising space utilisation in ORs with space constraints.

Benefits

Improved workflow efficiency

The floor-free Atlan workstations can reduce the time required for set-up and cleaning to improve room turnaround and streamline daily workflows.

Ceiling version:

- Easy positioning (motion control) and enhanced mobility to improve patient access.
- Built-in magnetic and optoelectronic coupling control to enable safe, secure and time-saving docking and undocking of the anaesthesia machine to the ceiling supply unit.
- Short connections between the anaesthesia machine and gas terminal units, power sockets and data terminals located on the media panels of the ceiling supply unit, to reduce risk of unintended disconnections of hoses and cables.
- Equipped with service trolley to facilitate safe and easy transport in case the machine must be maintained or exchanged.
- Additional mounting possibilities directly on the ceiling supply unit for workstation components (e.g. patient monitor, IT hardware, IV pumps, etc.) and accessories, to enable quick and easy exchange of the anaesthesia machine only, without the need to dismantle other components and accessories.

Wall version:

- Lifting mechanism for vertical positioning of the anaesthesia machine enabling precise working height adjustment with minimal physical effort.

Improved clinical workplace design and ergonomics

The floor-free Atlan anaesthesia workstations can improve the physical conditions of e.g. working environment, thereby contributing to performance, comfort, and job satisfaction.

Ceiling version:

- Broad range of horizontal and vertical positioning options to optimise space utilisation.
- Integrated cable ducts and cable channels for a well-organised, clutter-free working area, so as to reduce the risk of tripping hazards.
- Integrated grip handle, EasyGlide™ bearings and pneumatic brakes to enable precise and intuitive adjustment and placement of the working position (radius and height) for different users and surgical requirements.
- High carrying capacity of the ceiling pendant to accommodate various workstation components (e.g. with patient monitor, IT hardware, IV pumps, shelves, drawers, extension arms, and other components), increasing configurability of various approved workplace combinations to keep everything within reach.

Wall version:

- Optimises space utilisation in small ORs through its minimal footprint and depth.

Improved infection prevention and control

Recognising that medical devices can be part of the infection chain in hospitals, floor-free Atlan anaesthesia workstations are developed with infection prevention in mind to support hygiene standards within the OR.

- Designed with rounded edges, flat surfaces, and integrated cable management to comply with infection prevention standards.
- Channelling cables and wires via cable ducts and channels reduces potential contamination spots.

Benefits

Sustainability

Our technology protects people and the environment. Protecting the environment is therefore not just our duty, but also an opportunity for us to gain a competitive advantage through environmentally friendly products and by conserving resources. For this reason, we make sure that our day-to-day activities, in all areas of the company, protect the environment rather than damage it. We are committed to the United Nations' climate action goals. We handle resources responsibly and work systematically to reduce our carbon footprint. What is more, Dräger products and services help our customers all over the world live up to their responsibility to protect people and environment.

Learn more: https://www.draeger.com/en_uk/About-Draeger/Sustainability

System Components

D-15330-2024



Dräger Vapor 2000 and D-Vapor

As all other Dräger vaporisers, the Dräger Vapor 2000 series and D-Vapor deliver remarkable performance when it comes to precise agent delivery, safety, robustness, quality, and durability which may result in improved workflow efficiency, staff satisfaction, and clinical outcome.

D-30739-2017



Infinity Acute Care System

Transform your clinical workflow with Infinity Acute Care System. Its multiparameter monitor integrates with its networked medical-grade workstation, giving you real-time vital signs, access to clinical hospital systems and data management applications for a comprehensive range of patient information and powerful analysis tools at the point-of-care.

System Components


D-8549-2023



Vista 300

The Dräger Vista 300 is an essential patient monitor for adults, paediatrics, and neonates. It is part of a scalable patient monitoring system that fits securely into your hospital's network to give you all-time data access and optimise your workflows. Vista 300 facilitates enhanced workplace integration in your acute care areas with a comprehensive set of parameters and modules, helping you to improve patient care and outcomes as well as lower costs.

D-45302-2021

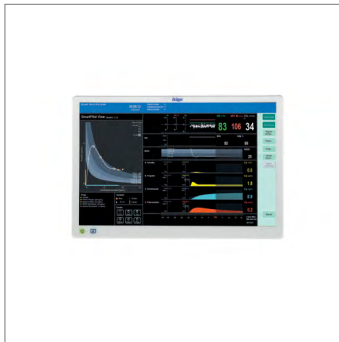

Dräger SmartPilot View

The software visualises the complex synergies of anaesthetic drugs and predicts their effects based on pharmacodynamic modelling for both the current status and the prospective course of general anaesthesia. SmartPilot View turns abstract device data into comprehensive visual information to support the decision on more precise and patient-optimised titration of anaesthetic drugs.

D-26017-2020



Dräger Ambia

Adapt acute care workstations ideally to your specific needs with our Ambia ceiling-mounted supply unit. With its numerous accessories and wide range of variation options, the Ambia ensures you maximum flexibility at the workplace. This not only helps to improve hospital workflows, but also increases the well-being of staff and patients alike.

Accessories

D-14586-2009



Infinity ID-accessories

Each and every Infinity ID-accessory has been designed to offer additional functionality, which can help you simplify routine tasks, streamline workflow and increase safety levels.

D-14348-2017



WaterLock 2

Protects your patient, protects your gas measurement systems. Designed to give you reliable gas measurements, the Dräger WaterLock 2 helps you to effectively filter humid and contaminated exhaled air thanks to our advanced membrane technology. Keeping your patients and investments safe from water, bacteria, and potential viruses.

D-6414-2018



Drägersorb 800+ – Soda Lime

Click and connect with 100% reliability. As one of the leading manufacturers of anaesthesia equipment, we believe in leading the way to producing high-quality soda lime that ensures your patients' and staff's safety to the highest degree. Drägersorb is more than just a formula, it is absorption efficiency — you can trust.

Related Products

D-12642-2018



Dräger Atlan A300/A300 XL

The new platform offers flexibility for most spatial conditions. The high precision piston ventilator supports lung protective ventilation measures and a comprehensive set of parameters assist decision-making support. The Atlan A300/XL can be networked to communicate securely with other networked devices to share data and information that can help to increase efficiency and reduce errors in anaesthesia.

D-2308-2022



Dräger Atlan A350/A350 XL

The new platform offers flexibility for most spatial conditions. The high precision piston ventilator supports lung protective ventilation measures and a comprehensive set of parameters assist decision-making support. The Atlan A350/XL can be networked to communicate securely with other networked devices to share data and information that can help to increase efficiency and reduce errors in anaesthesia.

D-316-2023



Dräger Perseus A500

The ceiling-mounted version of the Dräger Perseus A500 combines the award-winning, ergonomic design of the Perseus A500 with the mobility of a ceiling supply unit (CSU). In combination with Dräger Ambia, this floor-independent anaesthesia device offers a completely flexible and efficient workplace.

Technical Data

Dimensions of the ceiling and wall version

A tolerance of ± 10 mm (0.39 in) applies to all specified dimensions.

Dimensions of the ceiling-wall version (may deviate with accessory equipment)

Width	745 mm (29.3 in)
Height	1213 mm (47.8 in)
Depth	635 mm (25.0 in)

Dimensions of the service cart

Width	615 mm (24.2 in)
Height	576 mm (22.7 in)
Depth	625 mm (24.6 in)

Dimensions of the ceiling version with integrated castors (may deviate with accessory equipment)

Width	745 mm (29.3 in)
Height	1282 mm (50.5 in)
Depth	703 mm (27.7 in)

Weight

Nominal configuration consisting of electronically controlled gas mixer, plug-in connector for 1 vaporiser including vaporiser, breathing system, CLIC adapter and CLIC absorber, breathing hoses, central supply hoses (5 m (16.4 ft)), scavenging hose (5 m (16.4 ft))	approx. 120 kg (265 lbs)
Permissible total weight	300 kg (661 lbs)

Docking height to the ceiling supply unit

Ceiling-wall version with service cart Note: Docking height identical to Dräger Perseus A500	1021 mm (40.2 in)
Ceiling version with integrated castors Note: Docking height identical to Dräger Primus / Primus IE	870 mm (34.3 in)

Not all products, features, or services are for sale in all countries. Trademarks mentioned herein are the property of its respective owner. Trademarks may be owned by Drägerwerk AG & Co. KGaA (Dräger) or its affiliates in certain countries and not necessarily in the country in which this material is released. Visit www.draeger.com/trademarks for the current status of Dräger's trademarks.

Corporate Headquarters

Drägerwerk AG & Co. KGaA

Moislinger Allee 53–55
23558 Lübeck, Germany
■ www.draeger.com

Manufacturer

Drägerwerk AG & Co. KGaA
Moislinger Allee 53–55
23558 Lübeck, Germany

Region Europe

Drägerwerk AG & Co. KGaA

Moislinger Allee 53–55
23558 Lübeck, Germany
☎ +49 451 882 0
☎ +49 451 882 2080
✉ info@draeger.com

Region Middle East, Africa

Drägerwerk AG & Co. KGaA

Branch Office
P.O. Box 505108
Dubai, United Arab Emirates
☎ +971 4 4294 600
☎ +971 4 4294 699
✉ contactuae@draeger.com

Region Asia Pacific

Draeger Singapore Pte. Ltd.

61 Science Park Road
The Galen #04-01
Singapore 117525
☎ +65 6872 9288
☎ +65 6259 0398

Region Central and South America

Dräger Indústria e Comércio Ltda.

Al. Pucuruí 51 - Tamboré
06460-100 - Barueri - São Paulo
☎ +55 (11) 4689-4900
✉ relacionamento@draeger.com



Locate your Regional Sales
Representative at:
www.draeger.com/contact