

# Smarter anaesthesia with SmartPilot View

Keeping the right balance is a challenge when mixing anaesthesia drugs.

The software estimates and visualises complex anaesthetic drugs interactions, based on pharmacodynamic modelling for the current and predicted course of anaesthesia. SmartPilot View supports informed decision-making that facilitates more precision and individualisation in anaesthetic drugs titration.


Information from Dräger anaesthesia machines and IV-pumps is automatically transferred to SmartPilot View. It is plug and play. Additionally, manual entry is possible, e.g. for manual bolus.




## One dimensional view\* - Anaesthesia level over time

White line - the time course of the resulting level of anaesthesia, both for the past and the future 20 minutes.

Dashed line - appears when either anaesthesia parameter is changed but not confirmed yet or the “what if” function is selected. Based on the current pump rates and vapor settings it illustrates what would happen if a rate change or bolus were applied actually.

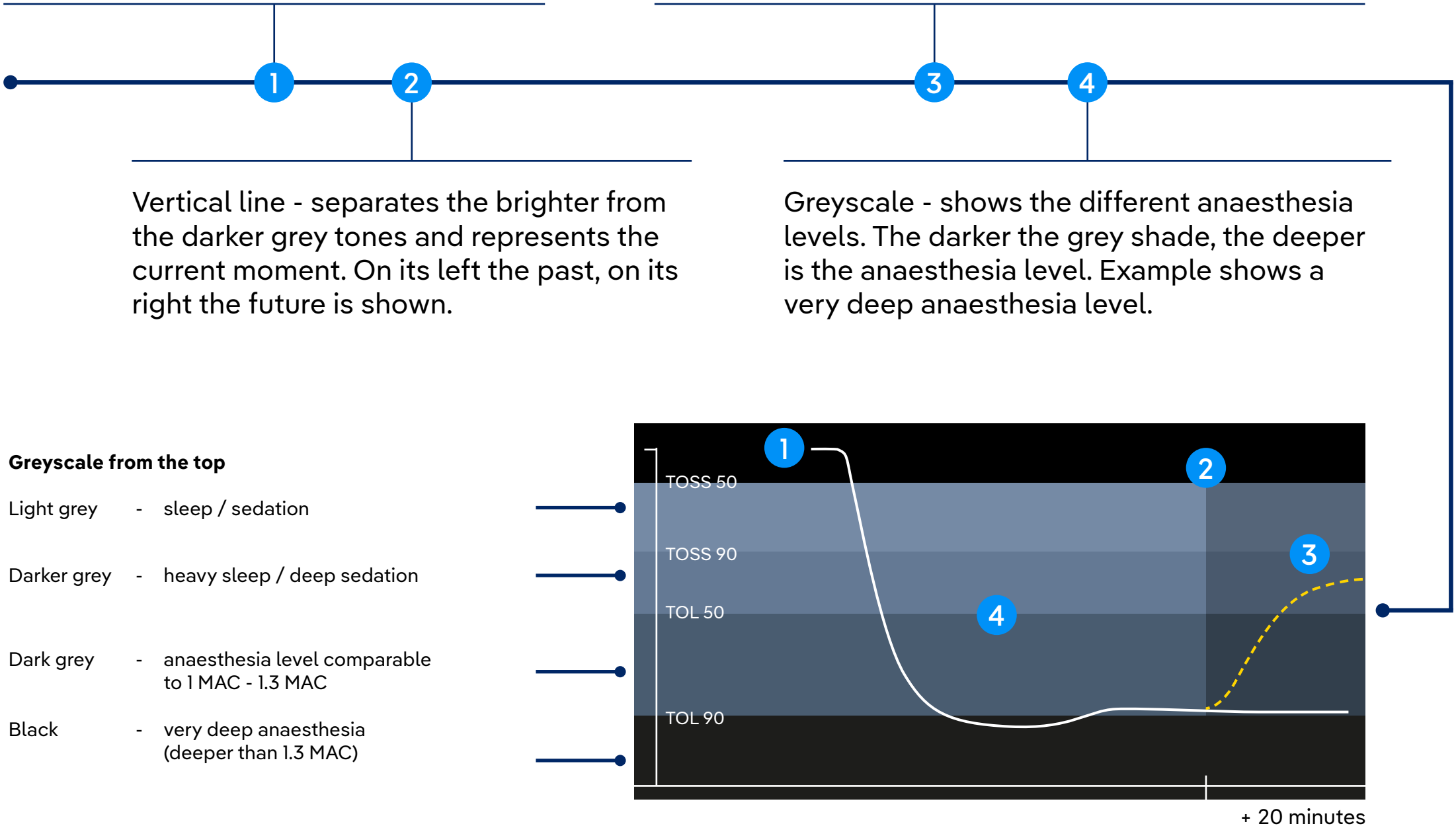

## Two-dimensional view\*\* - Anaesthesia level with visualised synergism of hypnotic and analgesic drugs

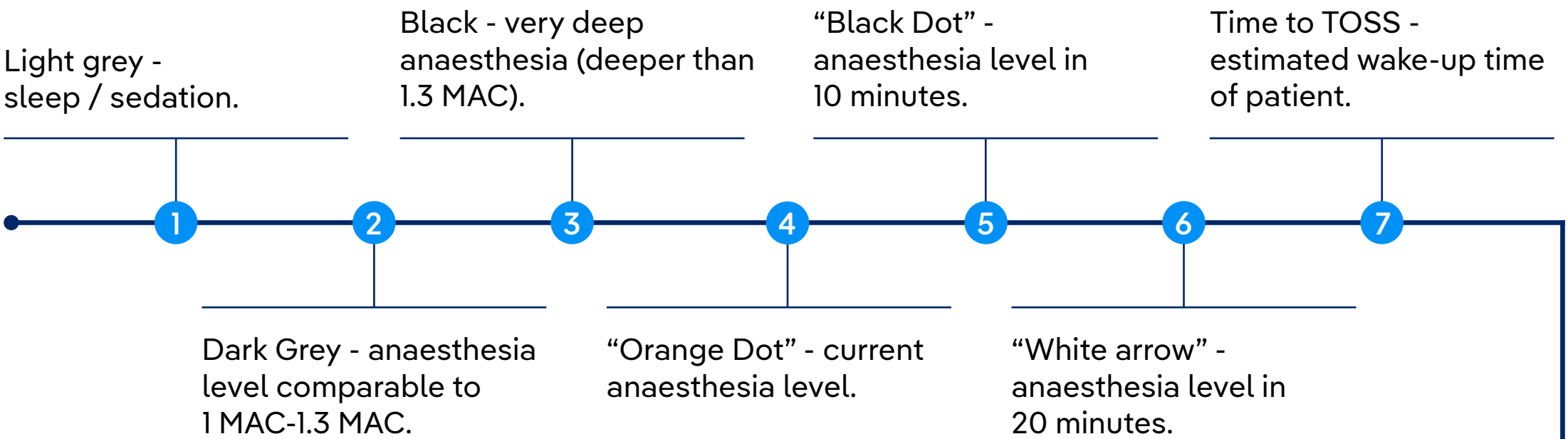


### Abbreviations

TOSS: Tolerance of Shake and Shout  
Example: TOSS 90 means 90% of patients do not react to shake and loud shouting

MAC: Minimum Alveolar Concentration. Used in mixed anaesthesia.  
Example: MAC 90 means 90% of patients do not react on surgical stimuli (comparable to 1.3 MAC); MAC 50 would be comparable to 1 MAC

TOL: Tolerance Of Laryngoscopy. Used in Total Intravenous Anaesthesia.  
Example: TOL 90 means 90% of patients tolerate a laryngoscopy

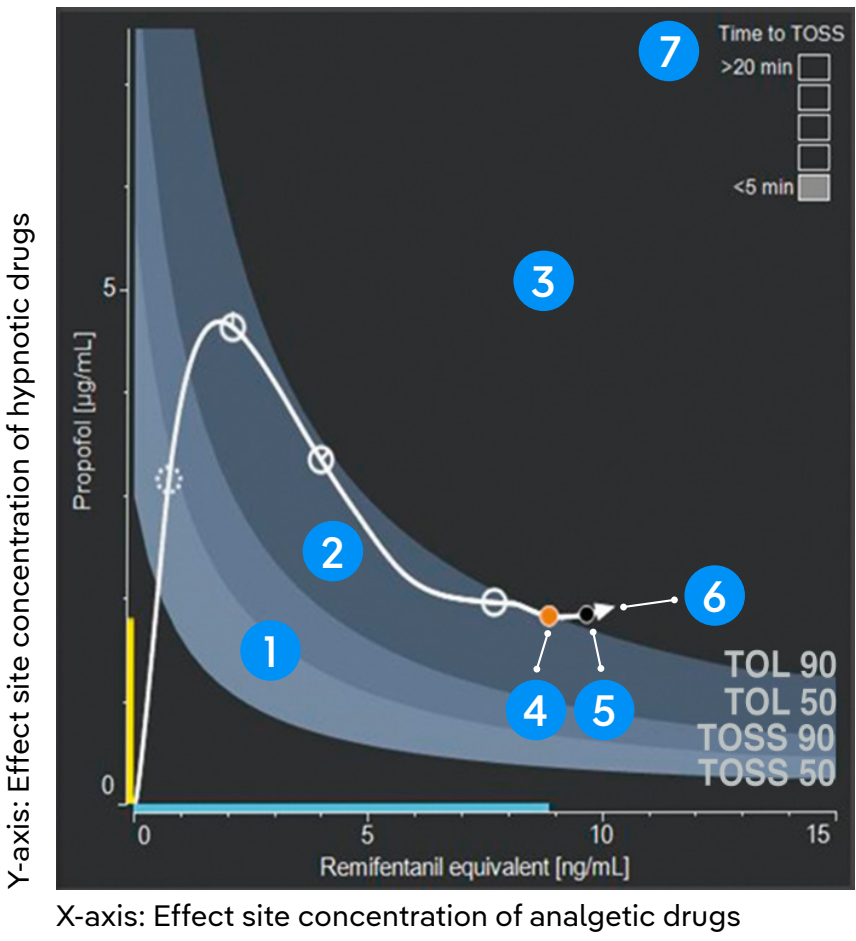

\* Noxious Stimulus Response Index  
\*\* 2D diagram for pharmacodynamics