



Dräger Perseus A500 Ceiling Anaesthesia Workstation

The ceiling-mounted version of the Dräger Perseus A500 combines the award-winning, ergonomic design of the Perseus A500 with the mobility of a ceiling supply unit (CSU). In combination with Dräger Ambia, this floor-independent anaesthesia device offers a completely flexible and efficient workplace.

Benefits

A well-established team

The maximum carrying capacity of Ambia allows for different approved workplace combinations with patient monitoring, IT systems and infusion pumps. The power supply units of Dräger monitoring and IT systems can be integrated into the device column of the ceiling-mounted version of the Perseus A500. Additional mounting options for components on the mounting column of the CSUs increase the configurability of the workplace.

Mobility on a high level

The floor-independent anaesthesia workstation makes best use of the available space. The double arm systems of Ambia have a large range so the workstation can be adjusted to the specific room requirements of the respective surgery. With EasyGlide™ bearings and pneumatic brakes the ceiling-mounted version of the Perseus A500 can be positioned exactly and with minimal effort.

Easy operation

One of the big advantages is the motion control. The workstation can easily be moved from the front side using the integrated grip handle. With this handle the anaesthesiologist can adjust the device position to new situations simultaneously to his patient care tasks. The attending physician can smoothly adjust the workstation height of the Perseus A500 to his or her body height. A safety sensor under the work surface interrupts the downward movement as soon as it registers an barrier, such as the legs of the sitting user.

Cable management

The Dräger cable management system helps in that cables can easily be laid in cable ducts. This way, all important connection and supply cables are protected from dust and contamination. Besides the Perseus A500 and the mounting arms for Monitoring and IT also the Ambia offers cable management solutions. The reduction of cable clutter not only reduces tripping hazards, but is also helpful from a hygienic perspective as the workplace design with its smooth, rounded surfaces helps cleaning and disinfecting the system.

Safety and transport

Built-in magnetic contacts ensure that the Perseus A500 remains mechanically coupled to the CSU if it is accidentally placed on a barrier. But if required, the Perseus A500 workstation can be disconnected from the CSU within a minute. The magnetic contacts of the system are deactivated via the key switch and the Perseus A500 can be placed on a special service trolley. As soon as the Perseus A500 is firmly fixed on this service trolley it can be transported comfortably.

System Components

D-7420-2011



Dräger Vapor 3000 and D-Vapor 3000

As all other Dräger vaporisers, the Dräger Vapor 3000 series and D-Vapor 3000 deliver same remarkable performance when it comes to precise agent delivery, safety, robustness, quality, and durability. In combination with the Perseus A500 anaesthesia workstation, they additionally can contribute to improved workflow efficiency, clinical outcome, and decision-making support.


D-30739-2017



Infinity Acute Care System

Transform your clinical workflow with Infinity Acute Care System. Its multiparameter monitor integrates with its networked medical-grade workstation, giving you real-time vital signs, access to clinical hospital systems and data management applications for a comprehensive range of patient information and powerful analysis tools at the point-of-care.

D-45302-2021

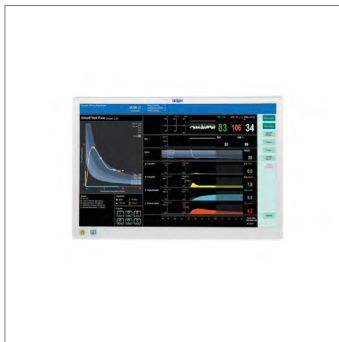

Dräger SmartPilot View

The software estimates and visualises complex interactions of anaesthetic drugs based on pharmacodynamic modelling for both the current and predicted course of anaesthesia. SmartPilot View supports informed decision-making that facilitates more precision and individualisation in anaesthetic drugs titration.

Accessories

D-14586-2009



Infinity ID-accessories

Each and every Infinity ID-accessory has been designed to offer additional functionality, which can help you simplify routine tasks, streamline workflow and increase safety levels.

D-14348-2017



WaterLock 2

Protects your patient, protects your gas measurement systems. Designed to give you reliable gas measurements, the Dräger WaterLock 2 helps you to effectively filter humid and contaminated exhaled air thanks to our advanced membrane technology. Keeping your patients and investments safe from water, bacteria, and potential viruses.

D-6414-2018



Drägersorb 800+ – Soda Lime

Click and connect with 100% reliability. As one of the leading manufacturers of anaesthesia equipment, we believe in leading the way to producing high-quality soda lime that ensures your patients' and staff's safety to the highest degree. Drägersorb is more than just a formula, it is absorption efficiency — you can trust.

D-19876-2015



FlexibilityTrolley

Dräger can support you already during OR planning, construction, or renovation periods. The Flexibility Trolley allows the Perseus A500 to be operated both as ceiling-mounted and as trolley-based workstation whenever required. This way, it offers the perfect bridging solution in hospitals where purchasing of anaesthesia equipment and OR infrastructure is planned consecutively. You can use the Perseus as trolley-based variant until the ceiling supply units are installed. It may even become a permanent solution.

Related Products

D-7118-2011



Dräger Perseus A500

Outstanding ventilator technology meets the latest approaches to ergonomics and system integration in one innovative anaesthesia machine, developed together with experts from all over the world to streamline your anaesthesia workflow.

D-45895-2021



Dräger Atlan A300/A350 Ceiling and Wall

The Atlan family's platform design offers floor-free anaesthesia workstations with ceiling and wall-mounted versions. They support organised clinical environments by enhancing workflow efficiency, improving patient access, simplifying cleaning procedures, and optimising space utilisation in ORs with space constraints.

D-9003-2016



Zeus Infinity Empowered

Lung-protective ventilation functionalities, cost-saving gas delivery technology, decision-support tools and comprehensive monitoring. Zeus IE integrates them all into one advanced anaesthesia workstation. It can, hence, improve clinical outcome, workflow efficiency, and cost savings throughout the anaesthesia process.

Technical Data

Perseus A500 ceiling mounted version

Dimensions (W x H x D)	1150 mm x 1045 mm x 790 mm 45.23 in x 41.14 in x 31.1 in
Weight Perseus A500 ceiling-mounted version	approx. 100 kg (220 lbs)
Nominal configuration	approx. 100 kg (220 lbs)
Maximum configuration*	200 kg (441 lbs)
* Refer to the instruction for use of the ceiling supply unit and observe the maximum load capacity.	

Mains power supply

Electrical mains connection	100 to 240 V ~ 50/60 Hz
Power consumption	typically 70 W, max. 400 W
Maximum power input	4 A

Storage areas

Big working surface, additional shelf arrangements (optional), lockable drawer (optional)

Not all products, features, or services are for sale in all countries. Trademarks mentioned herein are the property of its respective owner. Trademarks may be owned by Drägerwerk AG & Co. KGaA (Dräger) or its affiliates in certain countries and not necessarily in the country in which this material is released. Visit www.draeger.com/trademarks for the current status of Dräger's trademarks.

Corporate Headquarters

Drägerwerk AG & Co. KGaA
Moislinger Allee 53–55
23558 Lübeck, Germany
■ www.draeger.com

Manufacturer

Drägerwerk AG & Co. KGaA
Moislinger Allee 53–55
23542 Lübeck, Germany

Region Europe

Drägerwerk AG & Co. KGaA
Moislinger Allee 53–55
23558 Lübeck, Germany
☎ +49 451 882 0
☎ +49 451 882 2080
✉ info@draeger.com

Region Middle East, Africa

Drägerwerk AG & Co. KGaA
Branch Office
Dubai, United Arab Emirates
☎ +971 4 4294 600
☎ +971 4 4294 699
✉ contactuae@draeger.com

Region Asia Pacific

Draeger Singapore Pte. Ltd.
61 Science Park Road
The Galen #04-01
Singapore 117525
☎ +65 6872 9288
☎ +65 6259 0398

Region Central and South America

Dräger Indústria e Comércio Ltda.
Al. Pucuruí 51 - Tamboré
06460-100 - Barueri - São Paulo
☎ +55 (11) 4689-4900
✉ relacionamento@draeger.com



Locate your Regional Sales
Representative at:
www.draeger.com/contact